

Instruction Manual of Double-sided Medicine Trolley(English)



Warning-To reduce the risk of injury, user must read instructions manual carefully.



I. **Products Name:** Medicine Trolley

II. **Type No.:** MT750-01

III. **Application &Product Introduction**

Parameters of ABS Double-Sided Medicine Delivery Cart: (50 seats)

Dimensions: 750*480*930mm

1. Suitable for healthcare professionals to provide emergency care, dressing changes, and medication to patients.;

2. Main body

Primarily composed of aluminum, steel, and ABS engineering plastic structure, with four load-bearing columns and aluminum alloy columns.

3. ABS double-layer bottom injection molding technology, integrated handrails on both sides of the tabletop with a recessed design to prevent items from slipping. Three-sided stainless steel guardrails, transparent soft glass tabletop; two ABS trash cans for easy sorting and storage of trash; one document box and one L-shaped IV stand. One storage basket and two sharps containers.

4. Front: Central locking, five-tiered medicine tray. The interior is divided into morning, noon, and evening medicine sections, accommodating 25 people. One large medicine tray is 24cm high. The drawers are equipped with 3x3mm dividers, allowing for flexible storage of various medicines and equipment.

Back: Central locking system, five-tiered medicine tray, divided into morning/noon/evening 25-person slots.

5. Four deluxe silent casters, two with brakes and two without brakes, usable in any condition, with braking function.

IV. **Detailed Structural Diagrams of Each Part**

Functional Description of Each Part:

1. **Tabletop** -- For placing medical equipment such as laptops and electrocardiographs, not exceeding 100KG.
2. **Guardrail** -- To prevent items placed in the top panel from falling out due to movement.
3. **Top Panel** -- For placing medicines, saline bottles, etc. (not exceeding 30KG).
4. **Caster** -- With brake function.
5. **Hand-tight screw** -- Infusion stand lifting switch, infusion stand locking.
6. **Drawer Lock** -- Controls the closing of all drawers. Prevents loss of items.
Single drawer not exceeding 15KG.
7. **Medicine Box** -- Places appropriate items as needed (Set the partition limit diameter).
Single drawer not exceeding 15KG.
8. **Back Medicine Box** -- Places appropriate items as needed (Set the partition limit diameter).
Single drawer not exceeding 15KG.
9. **Large Drawer** -- Places appropriate items as needed (Set the partition limit diameter).
Single drawer not exceeding 15KG.
10. **Hand-tight screw** -- to secure the work surface.
11. **Sharps box** -- for used needles.
12. **Trash can** -- for temporary storage of waste other than medical waste (each weighing no more than 3kg).

V. Assembly Instructions

1. After opening the packaging, place the product on a clean, flat surface.
2. Remove all accessories, some are in the bottom drawer



3. Install the 4 casters:

- 1) Gently tip the cart onto its side or turn it upside down (it is recommended to place a soft cloth underneath) to expose the caster mounting bases on the underside of the frame.
- 2) Engage the caster brakes, align the threaded stems with the mounting holes on the frame, and tighten them securely.
- 3) Carefully return the cart to an upright position on the floor.
- 4) Testing: Verify that the brakes lock the wheels and prevent the cart from sliding; then, release the brakes to ensure the wheels rotate freely.

4. Attach the side-mounted document box, sharps container, waste bin, etc.(as the following picture)



VI. Product Safety Precautions

1. Do not stand on the top panel or use the product as a platform, as this may cause the vehicle to tip over and cause injury.
2. Do not heat the top panel or any other components. Avoid direct contact with boiling water above 60°C.
3. This product is equipped with 4 casters, 2 of which have brakes. If you cannot push it, please check if the casters are in the braked position.
4. When pushing, be careful not to let the casters squeeze or run over the power cord. Unplug the power cord when not in use for a long time.
5. Avoid storing in high temperature, humid, and dusty environments. Temperature: -10°C ~ 40°C, relative humidity: not exceeding 80%RH, atmospheric pressure: 760hPa ~ 1060hPa. Store in a well-ventilated indoor environment free of corrosive gases.

Product Medicine Box: The medicine box can be configured according to requirements and can be divided into various compartments of different sizes.

VII. Product Cleaning and Maintenance

When cleaning, wring out a soft cloth soaked in a diluted neutral detergent solution and wipe the product. Then wring out a cloth soaked in clean water to remove any remaining detergent. Finally, dry the product with a clean cloth.


- ※ Do not use volatile substances (thinners, volatile agents, gasoline).
- ※ Do not rinse directly with water; be careful of electric shock.
- ※ This product is not heat-resistant. For thorough cleaning, please use a low-temperature sterilization method.

VIII. Store

Packaged products should be placed in a dry, well-ventilated room without corrosive materials or heat source.



HEBEI WEBIAN MEDICAL INSTRUMENT TRADING CO.,LTD
A1-2703 SHIDAI BUILDING,GUANG AN STREET.CHANG AN
DISTRICT.SHIJIAZHUANG CITY.HEBEI PROVINCE.CHINA

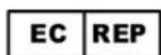

SUNGO Europe B.V
Fascinato Boulevard 522,unit 1.7.2909VA
Capelle aan del Ussel, The Netherlands



Quirumed S.L.U. CIF B97267405
C/ Corretger 117 Pq. Emp. Táctica, Paterna,
Valencia (Spain)

