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FOSHAN KAIYANG MEDICAL EQUIPMENT CO., LTD.
Dirección: 3 of No. 21. The 3rd Hongling Road, Shishan Town, Nanhai
District, Foshan City, Guangdong Province, China

SUNGO Europe B.V.
Fascinatio Boulevard 522, Unit 1.7, 2909VA Capelle aan den IJssel,
Países Bajos

Quirumed S.L.U. CIF B97267405
C/ Corretger 117 Pq. Emp. Táctica, 46980
Paterna, Valencia (Spain)

FOSHAN KAIYANG MEDICAL EQUIPMENT CO., LTD.

3 of No. 21. The 3rd Hongling Road, Shishan Town, Nanhai District, Foshan City, Guangdong Province, China

Medical Ward Screen

Instructions for use

Model: KY-P1/KY-P2/KY-P3

REF. 183-KY-P1S/183-KY-P2S/183-KY-P3S

Use:

The Medical Ward Screen is designed to separate a room or a part of a room or cubicle to provide temporary privacy for the occupant. It can be used around a bed, in front of a cubicle or in front of a window.

Composition:

The product is made of cloth (polyester fabric), steel round tube and wheels. Depending on the model, the screen ward consists of 1 screen, 2 or 3 screens.

Dimensions and weight:

- P1: 1800 × 1200 mm, N.W. – 4kg; G.W. – 5kg
- P2: 1800 × 1800 mm, N.W. – 6kg; G.W. – 7.5kg
- P3: 1800 × 2400 mm, N.W. – 8kg; G.W. – 9.5kg

Note: Each page is 1800 × 600 mm

Assembly:


1.Frame *1set

2. Square tube*2pcs


3.Wheel *6pcs

4.Small wrench *1pc

Install the tubes and attach the castors to the tubes:

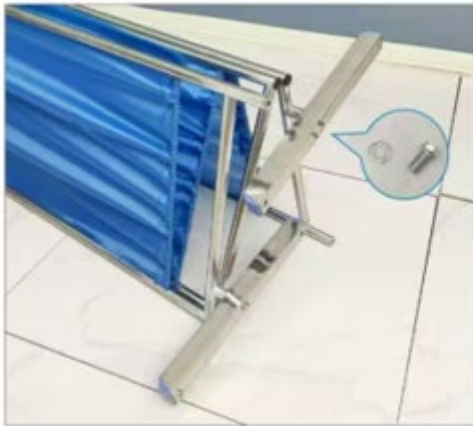Attach the castors to the round tubes on the left and right side of the frame, as shown below:

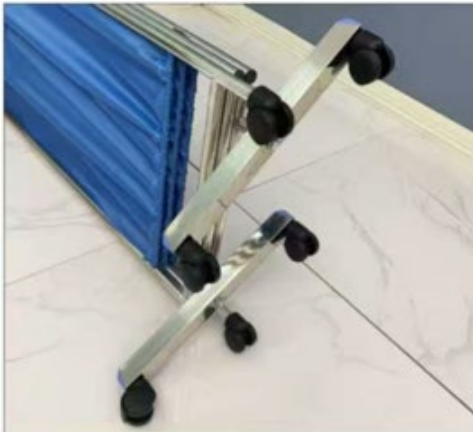Warning:

1. Assemble in the room of use.
2. Do not use this item if any parts are missing, damaged or worn.
3. Do not use this item unless all fixings are secured. Always use on an even surface. It is recommended that two people handle and assemble the item. Keep small parts out of reach of children.
4. Check if all screws or bolts are tightened and inspect regularly.
5. There are small components used in the construction of this product. These loose items should be kept away from children whilst assembling to avoid the danger of choking.

Cleaning methods:

Wipe with a neutral cleaner and wipe with a clean, soft cloth.

Storage:

1. Do not place the item on a wet or damp floor.
2. Periodically check all screws and fixings to ensure they are secure.

3. After unpacking the item, you may notice a smell due to the production process, however, this will disappear after a short period of time.

Transportation:


Avoid using the product under heavy pressure, direct sunlight, rain or snow.

Service life: 10 years

Disposal:


Please dispose of the product after use according to local regulations.

Symbols meaning:

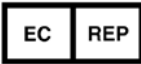Symbol	Introductions	Symbol	Introductions
	Batch Code		non-sterile
	Warnings and Precautions		Manufacture Date
	medical device		Consult instructions for use
	Use until year & month (Expiration date)		Name and Address of European Union Representative
	Manufacturer Name Address		CE Symbol
	Country of Manufacture		Model Number
	Importer		Distributor
	Unique Device Identification		

FOSHAN KAIYANG MEDICAL EQUIPMENT CO., LTD.
Address: 3 of No. 21. The 3rd Hongling Road, Shishan Town, Nanhai
District, Foshan City, Guangdong Province, China

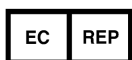SUNGO Europe B.V.
Fascinatio Boulevard 522, Unit 1.7, 2909VA Capelle aan den IJssel,
The Netherlands

Quirumed S.L.U. CIF B97267405
C/ Corretger 117 Pq. Emp. Táctica, 46980
Paterna, Valencia (Spain)

FOSHAN KAIYANG MEDICAL EQUIPMENT CO., LTD.

3 of No. 21. The 3rd Hongling Road, Shishan Town, Nanhai District, Foshan City, Guangdong Province, China

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SUNGO Europe B.V. Fascinatio Boulevard 522, Unit 1.7, 2909VA Capelle aan den IJssel, Niederlande