
User Manual

CE MD

Splint (EU)MDR2017/745

YXH-9 Basic UDI-DI:6974580870003DS

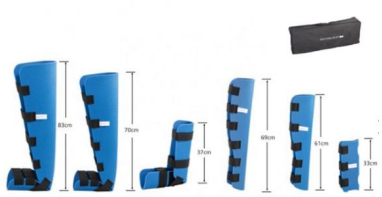I. Device description and specification

- This splint is made of EVA.
- Set of six size together.
- X-ray translucent and MRI compatible.

II. Intended Use

Splints are provided to people who need protection and support for painful, swollen or weak joints and their surrounding structures. Their designs make sure you position your Limbs（E.g. arm, finger, leg）correctly. It also provides support and stabilizer to the Limbs.

III. Instruction for Use

- Follow your doctor's instructions on how much weight you can put on your injured limb
- If the limb with the splint were not injured, wiggle them every now and then. This helps move the blood and fluids in the injured limb.
- If you have an adjustable splint that feels too tight, loosen it slightly.
- Keep up your muscle strength and tone as much as you can while protecting your injured limb. Your doctor may want you to tense and relax the muscles protected by the splint. Check with your doctor or your physiotherapist or occupational therapist for instructions.

IV. Cleaning and maintenance

Keep clean as routine. Wipe with 75% alcohol solution.

V. Storage

Ensure to avoid high temperature & moisture.

VI. Precaution and Warning

If your splint causes any of these problems, remove it and call your therapist right away: Pressure area such as sores or red marks that do not go away within one hour after removing the splint. / Increased swelling./ Excessive stiffness, pain or numbness./ Blisters.

If a splint gets wet, dry it with a hair dryer on the “cool” setting. Don’t use the warm or hot setting, because those settings can burn your skin.

Please dispose the product after use to comply with local regulation.

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Zhangjiagang Xiehe Medical Apparatus & Instruments Co., Ltd
EC REP **SUNGO Europe B.V.**