
Instruction Manual of
Three Function Hospital bed(English)

	Warning-To reduce the risk of injury, user must read instructions manual carefully.

I. **Products Name**: **Three function hospital bed**

II. **Type No**.: A02-6

III. **Structural composition**: (as picture)

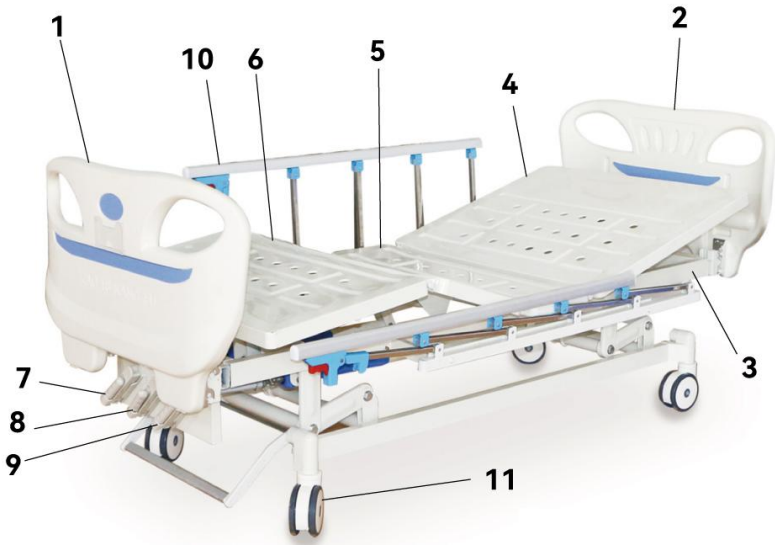





1) Bed Footboard
2) Bed Headboard
3) Bed-frame
4) Back panel
5) Welded bed panel
6) Leg panel
7) Crank for back lifting
8) Crank for overall lifting
9) Crank for leg lifting
10) Guardrails
11) Central brake Casters
12) Dining table
13) IV pole

IV. **Application:** It is suitable for patient nursing and recuperation, and facilitate the daily care to patient

1. The use of hospital beds should be monitored by professionals.
2. People who are taller than 2m and heavier than 200kg cannot use this bed.
3. This product should only be used by one person. Do not use two or more people at the same time.
4. The product has three functions: back lifting, leg lifting and overall lifting.

V. Installation:

1. Casters of bed

Brake the casters, then install the casters into the legs.

2. Bed headboard and footboard

The inner side of the headboard and the footboard is equipped with hanging inlay. The corresponding two metal mounting columns of the headboard and the footboard shall be pressed with vertical downward force to

embed the metal mounting columns into the inverted embedding groove, and locked with the hook of headboard and footboard.

3. Guardrails

Install the guardrail, fix the screws through the holes of guardrails and bed frame, fasten with nuts.

VI. **How to use：**

This hospital bed is equipped with three cranks, the functions are: back lifting, overall lifting, leg lifting.

1、Back rest lifting: Turn the crank clockwise, the back panel lift,
Turn the crank counterclockwise, the back panel down.

2、Overall lifting: Turn the crank clockwise, the overall lift,
Turn the crank counterclockwise, the overall down.

3、Leg rest lifting: Turn the crank clockwise, the leg panel lift,
Turn the crank counterclockwise, the leg panel down.

VII. **Attention**

1、 Check that the headboard and footboard were fasten tightly with bed frame.
2、 The safe working load is 120kg, the max load weight is 250kgs.
3、 After installing the hospital bed, place it on the ground and check whether the bed body shakes.
4、 The drive link should be lubricated regularly.
5、 Check the casters regularly. If they are not tight, please fasten them again.
6、 When operating the functions of back lifting, leg lifting and overall lifting, don't place the limb between the gap of bed frame and the bed panel or the guardrail, to avoid limb damage.

VIII. Transportation

The packaged products can be transported by general ways of transport. During the transportation, please pay attention to preventing sunshine, rain and snow. Avoid transportation with toxic, harmful or corrosive substances.

IX. **Store**

Packaged products should be placed in a dry, well-ventilated room without corrosive materials or heat source.

HEBEI WEBIAN MEDICAL INSTRUMENT TRADING CO.,LTD
A1-2703 SHIDAI BUILDING,GUANG AN STREET.CHANG AN DISTRICT.SHIJIAZHUANG CITY.HEBEI PROVINCE.CHINA

EC REP
SUNGO Europe B.V
Fascinato Boulevard 522,unit 1.7.2909VA
Capelle aan del Ussel, The Netherlands

Quirumed S.L.U. ClF B97267405
C/ Corretger 117 Pq. Emp. Táctica, Paterna,
Valencia (Spain)

CE MD